


EASY VALVE


Quick Start Manual




1. Switch heat on
(takes 3-5 min.)

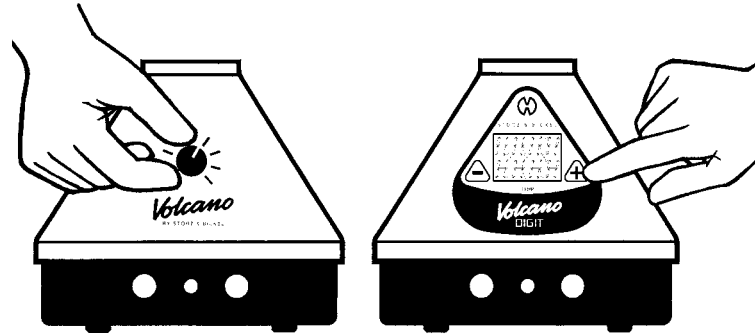


2. Choose temperature

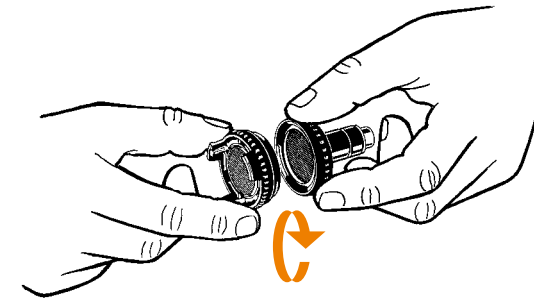


3. Remove Filling Chamber Cap

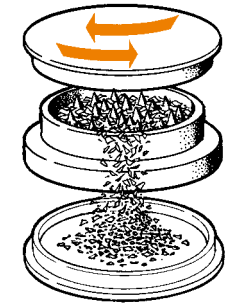


4. Grind substance

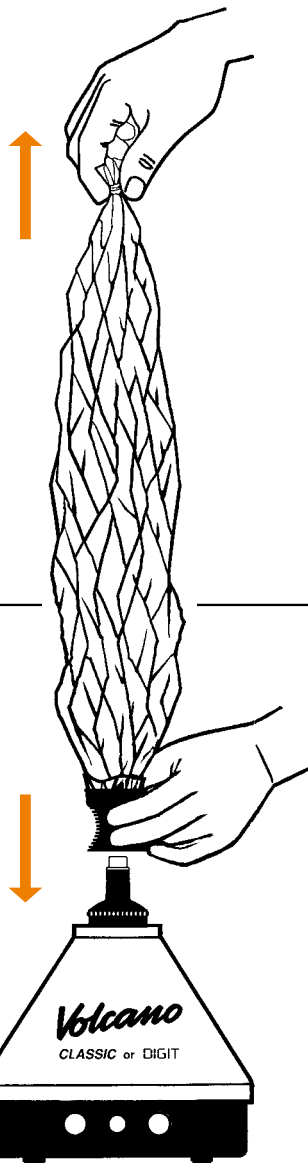

8. Before putting Valve Balloon on Filling Chamber, pull it tight

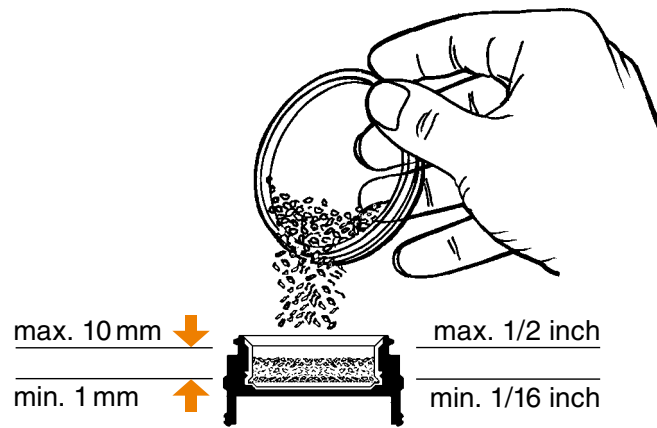


5. Fill Filling Chamber

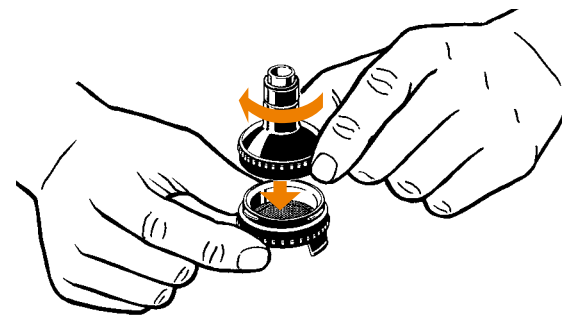

6. Screw on Filling Chamber Cap




7. Preheat substance if control lamp is off

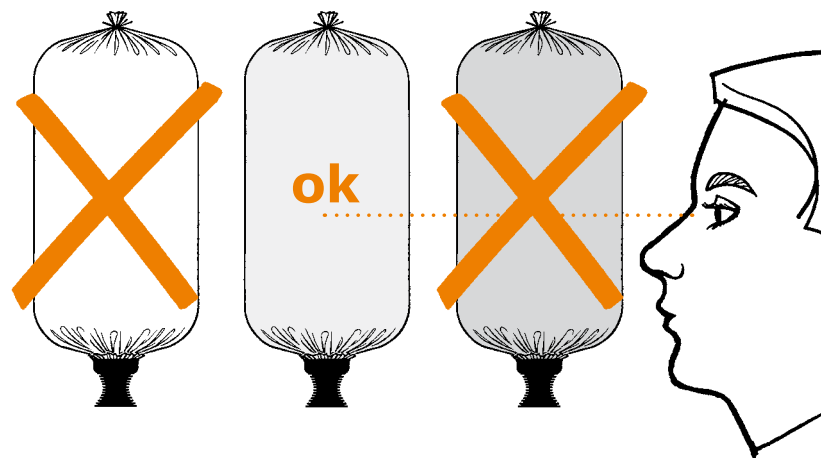

9. Watch vapor density




10. Remove full Balloon Unit

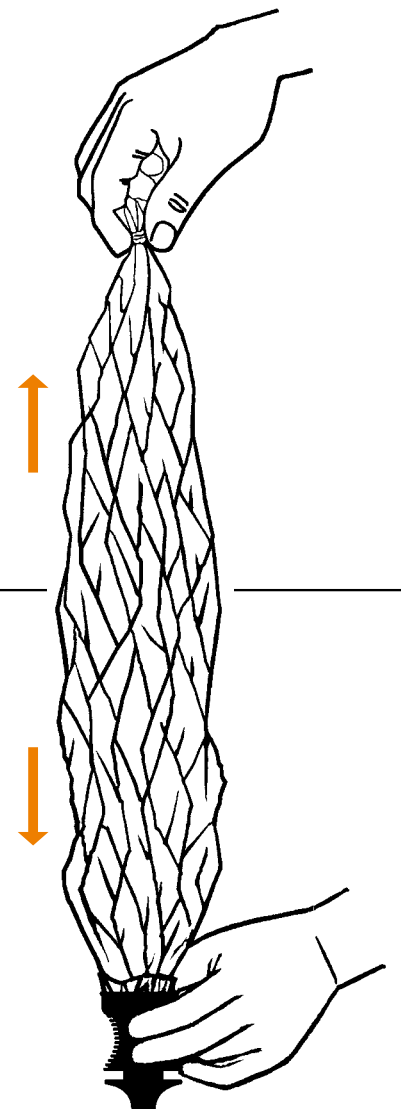


14. Empty the Valve Balloon completely by pulling it tight



11. Beware of hot surface!
Never leave Filling Chamber on VOLCANO

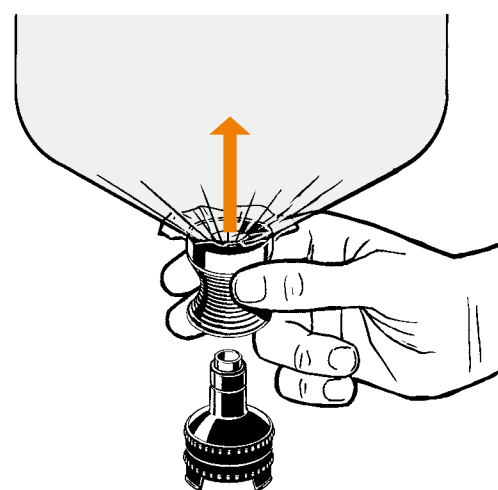

12. Detach Valve Balloon from Filling Chamber

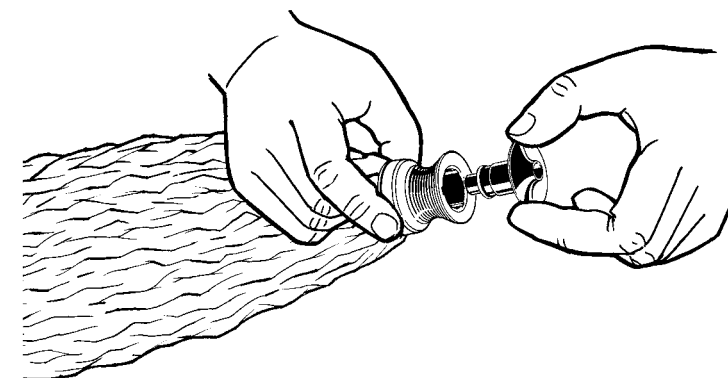

13. Attach Mouthpiece to Valve Balloon

The Instructions for Use of the VOLCANO Vaporization System are mandatory.